


Kromasil®

- Designed for ease of scale-up to preparative size columns
- High loadability
- Ultra-high purity spherical silica
- Superior durability
- ES Industries manufactures extremely high quality columns packed with Kromasil

Kromasil is manufactured by Eka Nobel, Sweden. It has been designed for both analytical and preparative scale liquid chromatography applications and is available in six particle sizes ranging from 3.5µm to 16µm. A choice of particle size enables the chromatographer to optimize a separation in the most economic manner. Kromasil's spherical uniform particles have a narrow pore size distribution and a unique surface chemistry. They give high chromatographic performance in both large and small scale applications. ES Industries has an extensive amount of experience in providing the best quality Kromasil packed columns in the industry.

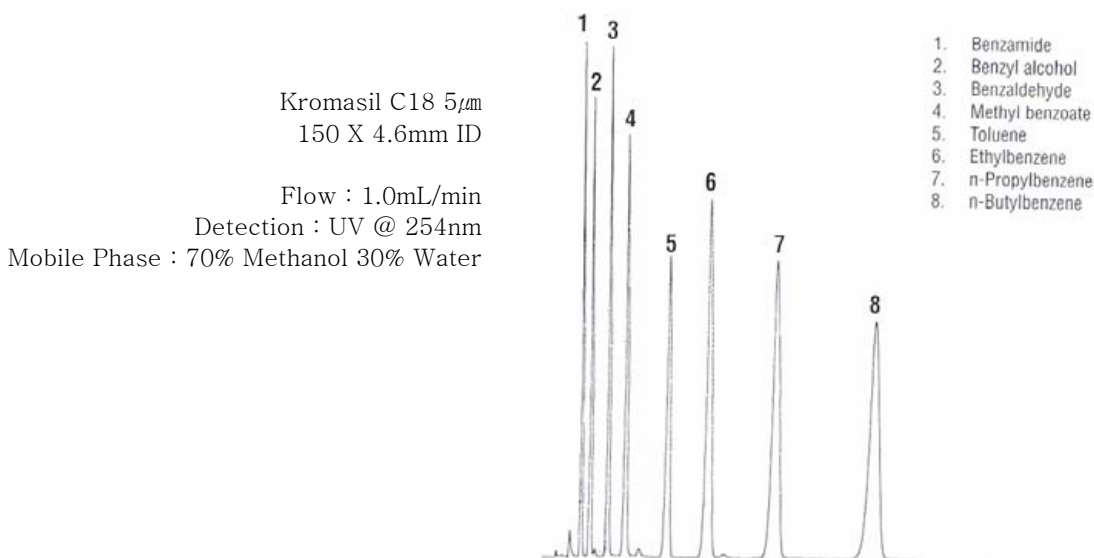

% Carbon Loading of Bonded Kromasil Phases

Material	Carbon %
Kromasil 100-C1	4.7
Kromasil 100-C4	8.0
Kromasil 100-C8	12.0
Kromasil 100-C18	19.0
Kromasil 100-NH2	1.5% Nitrogen
Kromasil 100-Silica	0.0